

Getting Started with NutriCurvist 1.0

! Note: Place the folder NutriCurvist_ver_01 directly on your desktop. Do not rename any files or folders, and do not move the folder into another location — this ensures the application runs correctly.

Data

NutriCurvist 1.0 uses a **built-in data entry table**

Use the **"Data Entry"** panel on the right side of the main window.

- Manually type or **paste** your data into the table.
- Each row represents one observation (replicate)

You can manage the table by right-clicking to paste data, add or remove rows, or use **Ctrl+V** to paste directly.

Tips:

For the analysis to run correctly, **every cell must have a valid entry**:

- No empty cells
- No missing or invalid numbers
- No text in numeric columns (X and Y)
- The first column (Sample_name) must contain at least one character, such as A1, A2, Sample1, or Sample2

If cells are empty or invalid, the app will either:

- Show an error message, or
- Return inaccurate or failed calculations.

Best Practices

- Use **at least 5 levels of nutrient and 4 replicates per group** for better modeling accuracy
 - Use the **"Clear Data"** button to reset the table
 - Always check your table for complete, valid entries before clicking **"Run Analysis"**
-

User guide for NutriCurvist

When starting **NutriCurvist.exe**, a window (Figure 1) will appear where the user is required to enter several values before running the model. These values are necessary for customizing the neural network and data processing pipeline for your specific dataset.

NutriCurvist
Version 1.0 (2025) | Machine learning powered tool for nutrient-response curve modeling

| Configuration

Number of bootstraps: *i*

Nutrient name: *i*

Response name: *i*

Unit of nutrient: *i*

Unit of response: *i*

Number of augmentations: *i*

Nutrient added noise (%): *i*

Response added noise (%): *i*

Custom nutrient level: *i*

| Data Entry

	Sample_name	X	Y
1	Sample 1	1.7400	13.7000
2	Sample 1	1.7400	25.7000
3	Sample 1	1.7400	29.8000
4	Sample 1	2.2400	36.2000
5	Sample 1	2.2400	29.3000
6	Sample 1	2.2400	34.1000
7	Sample 1	2.7400	52.0000
8	Sample 1	2.7400	56.5000
9	Sample 1	2.7400	56.4000
10	Sample 1	3.2400	81.5000
11	Sample 1	3.2400	85.8000
12	Sample 1	3.2400	80.1000
13	Sample 1	3.7400	87.9000
14	Sample 1	3.7400	76.6000

>> Run Analysis Help Contact Clear Data

University of Hohenheim | Institute of Animal Science

Figure 1. User interface of the NutriCurvist

Below is a description of each input field, what it means, what to enter, and how it is validated.

1. Number of Bootstraps
 - **What to enter:** A positive integer ≥ 10 (e.g., 100, 200)
 - **Meaning:** Determines how many random resampled datasets are created to ensure robustness and uncertainty estimation.
 - **Validation:** Must be a whole number (integer); 10 or higher
2. Nutrient name
 - **What to enter:** A text string (e.g., "Lysine", "Crude Protein")
 - **Meaning:** The name of the nutrient input variable; used as the X-axis label in plots.
 - **Validation:** Cannot be left empty

3. Response name
 - **What to enter:** A text string (e.g., "Weight Gain", "Egg Output")
 - **Meaning:** Name of the response variable; used as the Y-axis label in plots.
 - **Validation:** Cannot be left empty
 4. Unit of nutrient
 - **What to enter:** A text string (e.g., "%", "g/kg")
 - **Meaning:** The unit of measurement for the nutrient; shown in axis labels.
 - **Validation:** Cannot be left empty
 5. Unit of response
 - **What to enter:** A text string (e.g., "g/day", "mg/egg")
 - **Meaning:** Unit for the response variable.
 - **Validation:** Cannot be left empty
 6. Number of augmentations
 - **What to enter:** A non-negative integer (e.g., 0, 10, 50)
 - **Meaning:** Number of extra synthetic data points to generate using noise-added copies of your original data. Useful when working with limited datasets to improve model robustness. Set to 0 if you have sufficient data.
 - **Validation:** Must be 0 or higher
 7. Nutrient added noise (%)
 - **What to enter:** A non-negative number (e.g., 0, 2, 5)
 - **Meaning:** Amount of random noise (as a % of each original nutrient value) to be added to the nutrient-values for augmentation. Only applies if augmentation > 0.
 - **Validation:** Must be a non-negative number. If Number of augmentations is 0, then this must also be 0
 8. Response added noise (%)
 - **What to enter:** A non-negative number (e.g., 0, 5, 10)
 - **Meaning:** Amount of random noise (as a % of each original response value) to be added to the response-values for augmentation. Only applies if augmentation > 0.
 - **Validation:** Must be a non-negative number. If Number of augmentations is 0, then this must also be 0
 9. Custom nutrient level
 - **What to enter:** A non-negative number (e.g., 12.5) or the text 'na' to skip
 - **Meaning:** Allows scenario testing at a specific nutrient value. When set, the model evaluates the response at this nutrient level.
 - **Validation:** Must be a non-negative number or the text 'na'. If 'na' is entered, the analysis at a custom nutrient level is skipped
-

Results

Once NutriCurvist finishes analyzing your data and modeling is complete, the software automatically generates a results window with the following items inside:

This window contains **4 sheets**:

1. Metrics

- Includes metrics (e.g., R^2 , RMSE, nutrient required to reach desirable response, etc.) from the curve fitting process for each sample including their 95% confidence intervals.

2. Equations

- Contains the final symbolic formulas (ANN-based) fitted for each dataset. These are useful for interpretation or further use.

3. Coefficients

- Includes the estimated neural network parameters (weights and biases) for each sample:
- Each parameter also includes its **95% confidence intervals**, based on bootstrapping.

4. Data & predictions

Figures: Plots and metrics comparison

→ A **save** option is included in all sheets and Figures.

Insight from Results

The basic model used in **NutriCurvist** is a simple artificial neural network with **one artificial neuron** and a **tanh activation function**, designed to fit smooth and meaningful in the context of nutrient–response analysis (Figure 1 and Table 1).

Model Equation: The fitted model takes the following mathematical form:

$$\text{Response} = A \tanh(c \text{ Nutrient} + b) + B \quad \text{Equation 1}$$

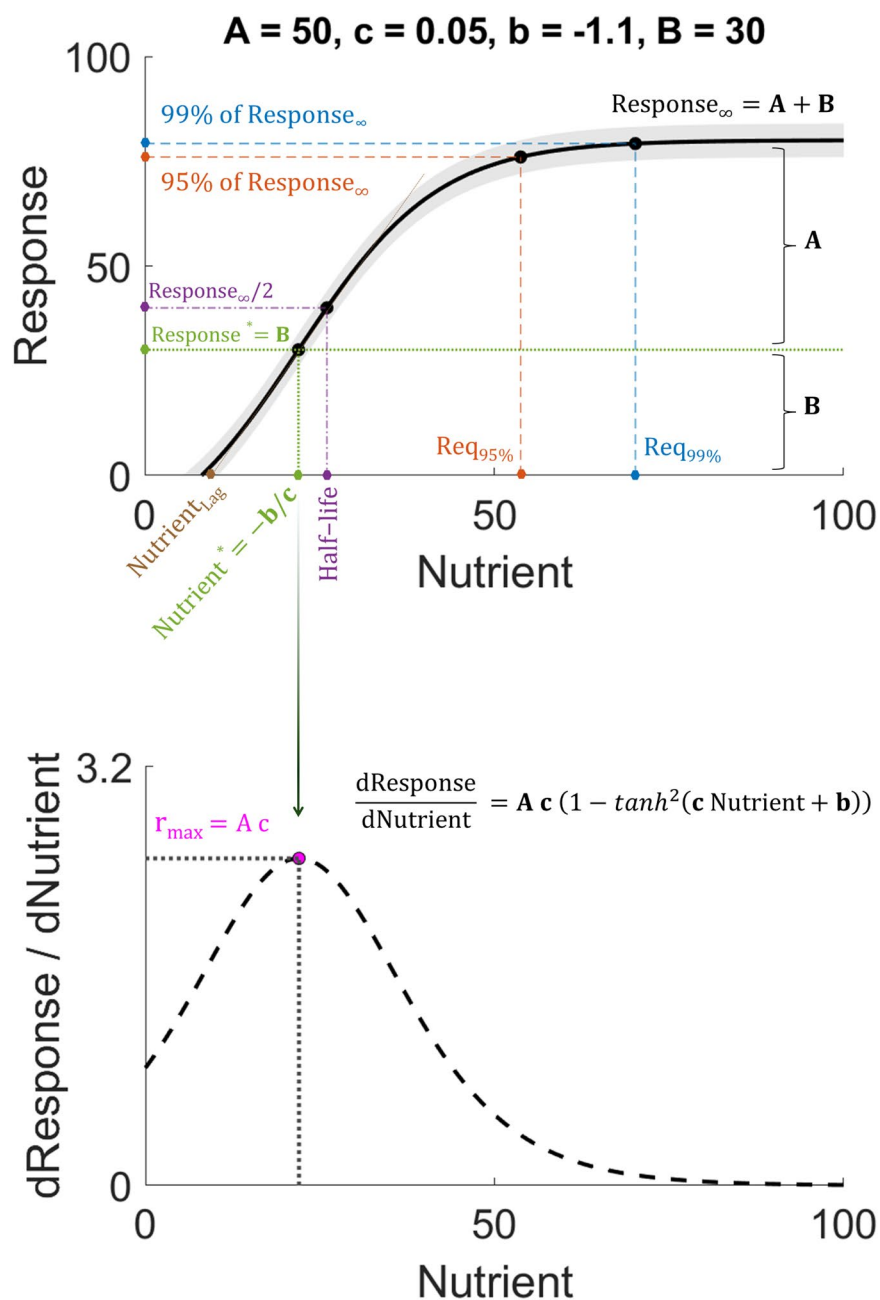


Figure 2: Schematic nutrient–response curve generated by the Equation 1, with parameters $A = 50$, $c = 0.05$, $b = -1.1$, and $B = 30$. For definitions see Table 1.

Table 1. Definition of NutriCurvist model, parameters, and derived nutritional metrics.

	Explanation	Equation
Model	Nutrient–response curve describing how response changes with nutrient	$\text{Response} = A \tanh(c \text{ Nutrient} + b) + B$
First derivative	Rate of change of the response with respect to nutrient	$\frac{d\text{Response}}{d\text{Nutrient}} = A c (1 - \tanh^2(c \text{ Nutrient} + b))$
Parameters		
A	Amplitude of the curve determines how much the response increases beyond the inflection point	
c	Steepness parameter: controls how quickly the response changes as nutrient level increases	
b	Horizontal translation: shifts the curve left or right along the nutrient axis	
B	Response at the inflection point; with $B \approx 0$, the curve rises symmetrically from the origin toward the plateau	
Nutritional metrics		
r_{\max}	Maximum slope of the curve; the steepest rate of increase in response per unit nutrient at the inflection point	$A c$
Nutrient*	Nutrient level at the inflection point, where response changes most rapidly	$-\frac{b}{c}$
Response*	Response value at the inflection point (i.e. Nutrient*)	B
Nutrient _{Lag}	Lowest nutrient value where a meaningful response begins (Response _{min} = model-predicted response at lowest nutrient)	$\text{Nutrient}^* - \frac{\text{Response}^* - \text{Response}_{\min}}{r_m}$
Half – life	Nutrient level at which 50% of the asymptotic response is reached	$\frac{\text{arctanh}\left(0.5 - 0.5\frac{B}{A}\right) - b}{c}$
Response _∞	Asymptotic response as nutrient increases indefinitely	$A + B$
Req _{95%}	Nutrient needed to achieve 95% of the Response _∞	$\frac{\text{arctanh}\left(0.95 - 0.05\frac{B}{A}\right) - b}{c}$
Req _{99%}	Nutrient needed to achieve 99% of the Response _∞	$\frac{\text{arctanh}\left(0.99 - 0.01\frac{B}{A}\right) - b}{c}$